


Vanderlift™ II Components

- A. Model Name Vanderlift™ II
- B. Four Point Floating Hanger Assembly The 4 point self adjusting sling attachment hanger assembly is standard on Vancare lifts. A two point and 4 point clip style is also available.
- C. Safety Latch Prevents the loop of the sling From accidentally coming off.
- D. Actuator Used to raise and lower the lift arm.
- E. Front Caster Dual 3" Caster is hospital grade.
- F. Rear Caster The heavy-duty 4" locking Caster is hospital grade.
- G. Locking Device Easy to lock/unlock with foot.
- H. Manual Base Spread Handle Used to open/close legs of the lift. Power leg spread option available.
- I. Control Box Houses electronics including battery gauge, emergency down button and emergency stop button.
- J. Emergency Stop Button Used to shut-off power to the lift. Normally in open (out) position.
- K. Battery Pack Can be charged on the lift or in a wall mounted charging station.
- L. Handle Assembly Used to maneuver the lift.
- M. Hand Control Used to raise and lower the lift as well as to spread the legs with the power leg spread base motor option.
- N. Manual Emergency Lowering Device Used to manually lower the lift in the event of power failure.
- O. Mast Joint Steel mast joint is more durable than plastic.



VANCARE®
We move body and mind

Vanderlift™ II
B 450, B 600

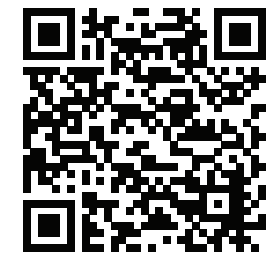

Scan the barcode above with a smart phone for more product information.

VANCARE, INC
1515 1st STREET
AURORA, NE 68818
TEL 800.694.4525
TEL 402.694.4525
FAX 402.694.3994
vancare.com



Made with Pride in the USA

You can trust Vancare

Trust Vanderlift™ II to provide safe transfer of residents/patients, regardless of range-of-motion limitations. The B450 model (shown on page 1) has a 450 lbs. capacity, while the B600 model (shown on page 2) serves residents/patients weighing up to 600 lbs.

Both models are proven to reduce back injuries common among nursing staff - helping to minimize lost work days, worker's compensation claims, staff turnover, and liability insurance premiums of service and use.



Vanderlift™ II B 450 and B 600 features

- Magnetized hand-held control secures anywhere on the lift.
- Floating four-point sling connection is easy to use and conforms to the body shapes for additional resident/patient comfort.
- Auxiliary hand-operated control provides means of lowering the resident/patient. (In case of power failure, the resident/patient can be safely lowered with an emergency manual system.)
- Removable, rechargeable 24 volt DC memory-free battery pack features a low-battery power indicator, and actuator.
- Optional built-in scale available.
- Features an optional electric expandable base for added versatility.
- Lifting height adapts from the floor to the highest rehab beds and tables.
- Optional Intellipac for keeping track of use and service.



The versatile Vanderlift™ Series Lift is the ultimate risk manager!

A versatile Vanderlift™ II Series Lift can easily adapt to most transferring situations, from lifting a resident/patient from the floor to raising him or her to the highest specialty rehab table. With models capable of handling patients/residents up to 450 or 600 lbs., even the largest and least ambulatory people can be transferred safely comfortably.

Patient Sling Features

Vancare's sling selection is world class. With several hundred models, sizes and material choices, Vancare can provide the perfect fit for your patients and your budget. Vancare also proudly produces top quality "Made in the USA" products. When your patient's safety is at stake, you can trust a Vancare sling.

